


TOYDAD 600 L

TWO CHANNEL UV DETECTOR
FOUR CHANNEL UV DETECTOR
SCANNING UV DETECTOR

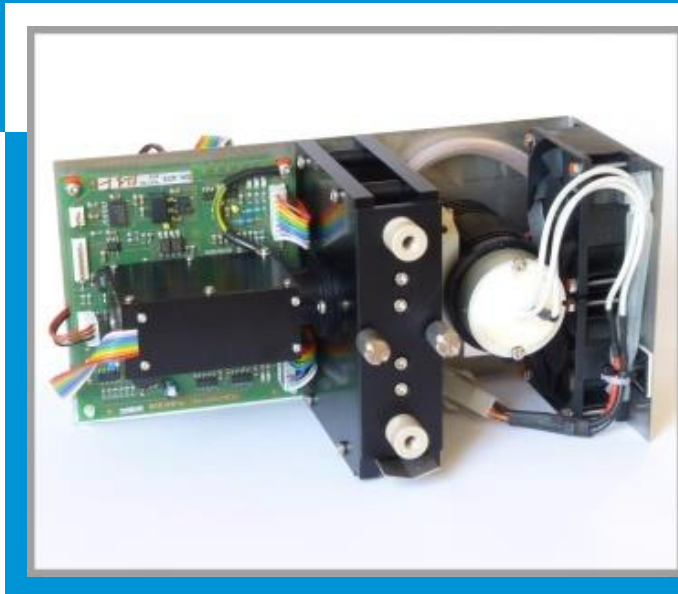

is **OEM built-in** UV detector suited for preparative and flash chromatography with diode array (DAD) in range 200 – 600 nm. It is offered in the following versions:

Two Channel version measures at **two wavelengths simultaneously** with possibility to see **current scan**.

Four Channel version measures at **four wavelengths simultaneously** with possibility to see **current scan**.

Scanning version measures at **four wavelengths simultaneously** or **sends scan with speed 5 Hz** which allows to **create 3D picture**.

More wavelengths or scan are used in liquid chromatography to **verify purity of analyzed samples** or in situations when some peaks absorb on different wavelengths.

Most important features of TOYDAD detectors:

- small dimensions
- **built-in lamp working hours counter**
- **easy replacement of flow cell**
- **measuring at more wavelengths simultaneously**
- **sophisticated diagnostic software**

SPECIFICATIONS

TECHNICAL PARAMETERS:

	TOYDAD 600 L Two Channel UV Detector	TOYDAD 600 L Four Channel UV Detector	TOYDAD 600 L Scanning UV Detector
Part Number	TOYL1100	TOYL1200	TOYL1300
Wavelength range	200 - 600 nm (256 elements on CCD)		
Number of channels (Signals)	2	4	4
Scan	Informative 200 – 600 nm	Informative 200 – 600 nm	200 - 600 nm, 5 Hz, step 2 nm
Typical spectral half-width	8 nm		
Accuracy of adjustment / Reproducibility	± 1 nm /±0.5 nm		
Noise level at test cell (254 nm, TC 0.75s)	± 5 x 10 ⁻⁵ AU		
Drift at test cell (254 nm after 1 h)	1 x 10 ⁻³ AU/hr		
Materials in contact with mobile phase	PTFE; fused silica, stainless steel		
Time constant (T63)	0.5 s, 0.75 s, 1.0 s, 2.0 s, 4.0 s, 8.0 s, 16.0 s, 0.2 s, 0.1 s		
Interface	RS232		
Power supply	24 -36 V DC		
Power input	48 W		
Dimensions (W x H x D)	220 x 120 x 110 mm (8.67 x 4.72 x 4.33 in)		
Weight	1.90 kg (4.19 lb)		

PREPARATIVE CELL PLCC 15 (SUPPLIED WITH UNIT)

volume/optical path (adjustable)	45 µl / 0.3 mm; 55 µl / 1.4 mm; 70 µl / 2.4 mm
Cell connecting	tubing with OD = 1/8", thread 1/4"-28
Maximal flow rate	500 ml/min